

Absorbest Super Foam


Absorbest Super Foam is indicated for use on exuding wounds. Absorbest Super Foam is a sterile, polyurethane foam dressing with the exudate management properties of a superabsorbent dressing. The wound dressing has a soft adhesive contact surface to facilitate dressing application.

Absorbest Super Foam is a thin, discrete wound dressing that can handle small to large amounts of wound fluid. The dressing locks in wound exudate more efficiently than conventional foam dressings and can hold thick wound fluid. The surrounding skin is prevented from maceration, while a good moisture balance is achieved in the wound area. The wound exudate is effectively transported away from the wound and distributed evenly into the dressing's core constructed by our DryMax Technology. Absorbest Super Foam absorbs wound fluid that contains bacteria and other harmful substances.


A soft adhesive layer covers the wound contact surface. This will enhance the dressing application and work as the caregiver's extra hand. The soft adhesive enables reposition and provides a gentle surface at dressing removal.

The blue hydrophobic backing reduces the risk of wound fluid strike-through. The Absorbest Super Foam is designed to work well together with compression therapy.


Application guide

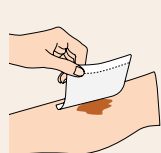

The dressing must be placed close to the wound, or placed in close combination with a suitable dressing.

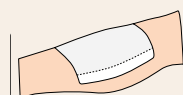

Secure the dressing with suitable bandage, film or tape.

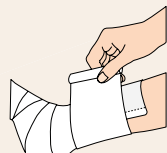

Absorbest Super Foam can be used in conjunction with compression therapy.

Item	Product size (cm)	REF	Pcs/box	Absorption capacity*
Absorbest Super Foam	7.5x7.5	F60041/10	10	60 ml
Absorbest Super Foam	10x10	F60043/10	10	130 ml
Absorbest Super Foam	15x15	F60044/10	10	350 ml
Absorbest Super Foam	10x20	F60045/10	10	350 ml
Absorbest Super Foam	20x20	F60046/10	10	640 ml

The maximum free swell absorption capacity* of Absorbest Super Foam is 1.8 ml per cm² of the product core.

*) According to international standard EN 13726



Manufactured by Absorbest AB Sweden
Distributed by Alliant Biotech, LLC.
www.alliantbiotech.com
888-307-1144



- **Cost effective**
- **User friendly**
- **Increased patient comfort**
- **Improved wound healing***

*compared to conventional wound dressings



DryMax® technology



Manufactured by Absorbest AB
Sweden
www.absorbest.com