

## DryMax Soft

DryMax Soft is a sterile, conformable, superabsorbent dressing for exuding wounds. Its properties provide an easy-to-use and reliable solution for both patient and clinician. DryMax Soft absorbs and retains exudate in the dressing. The exudate is absorbed vertically by the dressing and retained in the core, forming a gel, thus minimizing the risk of maceration and maintaining a moist and balanced wound environment.

The seamless cover simplifies positioning over the wound. The design allows the dressing to be applied with either side facing the wound. This feature also makes it possible to stack dressings on each other if extra absorption capacity is required.


The dressing remains dry and comfortable for the patient. Simultaneously, it contributes to a moist wound environment. The hydrophilic contact layer transports fluid effectively away from the wound and distributes it evenly into the dressing's core constructed by our DryMax Technology.

The ultra-thin and highly conformable features allow the dressing to be used on all body parts. DryMax Soft is designed to work well together with compression therapy.


### Application Guide



The dressing must be placed close to the wound, or placed in close combination with a suitable dressing.

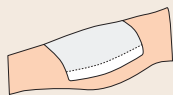

Secure the dressing with suitable bandage, film or tape.

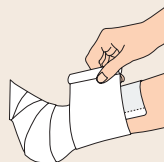

DryMax Soft can be used in conjunction with compression therapy.

| Item                     | Product size inches | REF       | Pcs/box | Absorption-capacity* |
|--------------------------|---------------------|-----------|---------|----------------------|
| DryMax Soft 7.5x7.5 cm   | 3" x 3"             | F60011/10 | 10      | 60 ml                |
| DryMax Soft 10x10 cm     | 4" x 4"             | F60013/10 | 10      | 130 ml               |
| DryMax Soft 10x15 cm     | 4" x 6"             | F60012/10 | 10      | 300 ml               |
| DryMax Soft 12.5x12.5 cm | 6" x 6"             | F60123/10 | 10      | 220 ml               |
| DryMax Soft 10x20 cm     | 4" x 8"             | F60015/10 | 10      | 330 ml               |
| DryMax Soft 15x20 cm     | 6" x 8"             | F60014/10 | 10      | 400 ml               |
| DryMax Soft 20x20 cm     | 8" x 8"             | F60016/10 | 10      | 600 ml               |
| DryMax Soft 20x25 cm     | 8" x 10"            | F60018/10 | 10      | 780 ml               |
| DryMax Soft 20x30 cm     | 8" x 12"            | F60017/10 | 10      | 960 ml               |
| DryMax Soft 37x56 cm     | 15" X 22"           | F60019/10 | 10      | 3000 ml              |

\*According to international standard EN 13726



- Cost effective
- User friendly
- Increased patient comfort
- Improved wound healing\*

\*Compared with conventional wound dressings



DryMax® technology



Distributed by Alliant Biotech, LLC.  
888-307-1144  
www.alliantbiotech.com



Manufactured by Absorbest AB  
Sweden  
www.absorbest.com